

Water bath

Shaking Water Bath

(진탕항온수조)


Feature



C-SKW1




Drain Valve

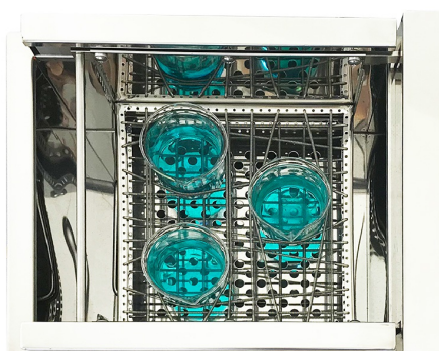

좌, 우 왕복

기본

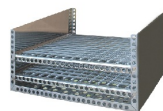

Spring Rack

시험관 크기에 상관없이 고정가능


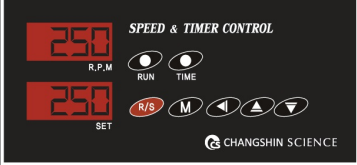
옵션




Flask Holder

시험관 크기별 Holder를 장착해 고정가능

Control Feature

	Temp Control <ul style="list-style-type: none"> - Temp range : Ambient ~ 100°C - Time range : 99 Hours 59 Minutes - Digital safety thermostal. - Auto tuning - Temp calibration - Digital P.I.D
	Speed Control <ul style="list-style-type: none"> - R.P.M : MAX 250 - 3.5cm Stroke - Time range : 99 Hours 59 Minutes

동 식물 세포의 진탕배양, 미생물 배양, 효소반응, 혈청검사, 임상병리, 각 종 효소측매의 반응 실험 등에 사용됩니다.

- ▶ 좌, 우 왕복 Shaking
- ▶ Max 250 R.P.M, 3.5 Stroke
- ▶ 소음이 적고 내구성이 우수한 D.C Motor
- ▶ Spring Rack 기본 장착 (Flask Holder 교체가능)
- ▶ Temp, Speed 개별제어 (Digital P.I.D)
- ▶  Arch Type Bath Cover 로 응축수가 실험물에 영향을 주지 않음.

Technical Data

Model		C-SKW1
Dimensions	Inner (W×D×H)(mm)	400×300×200
	Outer (W×D×H)(mm)	700×410×380
	Volume (L)	24
Controller		Digital Temp Controller & Digital Display / PT 100Ω
Temp	Range	Ambient ~ 100 °C
	Accuracy	0.3°C
Timer	Scale	99 hours 59 minutes
Electrical spec	Heater	2000W
	Power	220V, 9A, 1ø
Materials	Inner	Stainless steel SUS304
	Outer	Steel plate with powder coated
Safety device		Exclusive over temp. protector
Weight (Kg)		32